



GEMINI椎间融合系统

INTERBODY FUSION CAGE

LF-III型椎间融合器

LF-III型椎间融合器

INTERBODY FUSION CAGE




适应症：

- 退变性腰椎间盘疾病与腰椎不稳
- 腰椎间盘突出术后翻修术
- 退变性、峡部裂型脊柱滑脱 (I°-II°)

产品信息：

型号	规格	产品编号
LF-III	28×10×7	105020328007
LF-III	28×10×8	105020328008
LF-III	28×10×9	105020328009
LF-III	28×10×10	105020328010
LF-III	28×10×11	105020328011
LF-III	28×10×12	105020328012
LF-III	28×10×13	105020328013
LF-III	28×10×14	105020328014
LF-III	28×10×15	105020328015
LF-III	28×10×16	105020328016
LF-III	31×12×7	105020331207
LF-III	31×12×8	105020331208
LF-III	31×12×9	105020331209
LF-III	31×12×10	105020331210
LF-III	31×12×11	105020331211
LF-III	31×12×12	105020331212
LF-III	31×12×13	105020331213
LF-III	31×12×14	105020331214
LF-III	31×12×15	105020331215
LF-III	31×12×16	105020331216

产品特点：

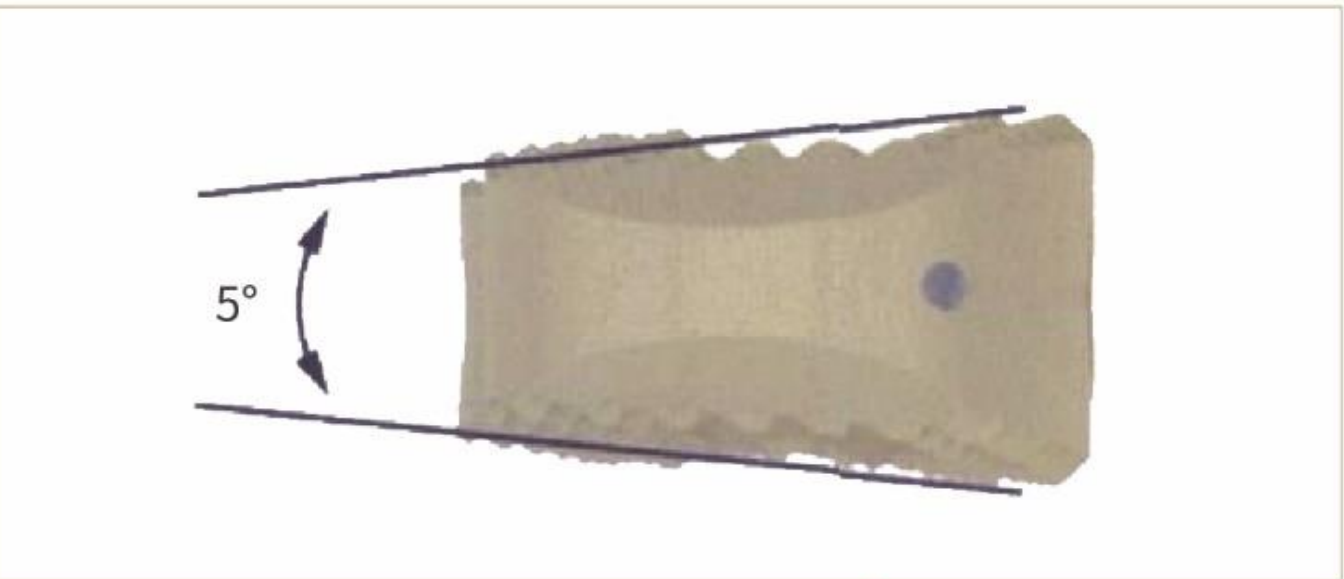

5°前凸角恢复生理前凸

锥形齿设计, 有效防滑, 提供初始稳定

钽标记, 射线下可见, 保证植入安全

超大植骨窗, 提高融合率

前端圆钝设计, 易于植入

PEEK材料, 生物相容性好, 弹性模量与皮质骨相近



浙江德康医疗器械有限公司
ZHEJIANG DECANS MEDICAL DEVICES CO., LTD

📍 浙江省嘉兴市秀洲区加创路321号上海交大(嘉兴)科技园7-8号楼
☎ 0573-82792526 📠 0573-82792536-8088 🌐 www.decansmd.com

